



# NONIN EQUANOX 8004CA ADVANCE CEREBRAL OXIMETER SENSOR PROVIDES VALID ASSESSEMENT OF TRUE TISSUE OXYGEN SATURATION



David B. MacLeod, FRCA, Keita Ikeda, Ph.D

Department of Anesthesiology, Duke University Medical Center, Durham, NC 27710, USA

## Introduction

Near infra-red spectroscopy (NIRS) regional oximeters estimate cerebral tissue oxygen saturation (rSO<sub>2</sub>) and is a composite value based on the relative proportions of arterial and venous blood within brain tissue. Nonin has developed a regional oximeter with dual emitter - detector sensor architecture to eliminate the influence of variations in non-cerebral tissue (eg:scalp). The aim of this study was to calibrate and validate a new four wavelength cerebral oximeter sensor.

## Methods


Adult ASA 1 volunteers were enrolled. An internal jugular venous bulb catheter and a radial artery catheter were placed. Cerebral sensors were placed bilaterally on the forehead and a pulse oximeter was placed on the ear.

NIRS-derived rSO<sub>2</sub> and pulse oximeter saturation (SpO<sub>2</sub>) were recorded continuously at 1Hz. Hypoxia was induced and managed via a dedicated facemask and breathing apparatus (Respiract, Thornhill Research, Toronto) and monitor continuous end-tidal CO<sub>2</sub> and O<sub>2</sub>.



Subjects underwent two standardized breath-down protocols in sequence:

- (1) decrease of SpO<sub>2</sub> in approximately 5% increments from 100 to 70%
- (2) decrease of SpO<sub>2</sub> in single step from 100 to 70%.

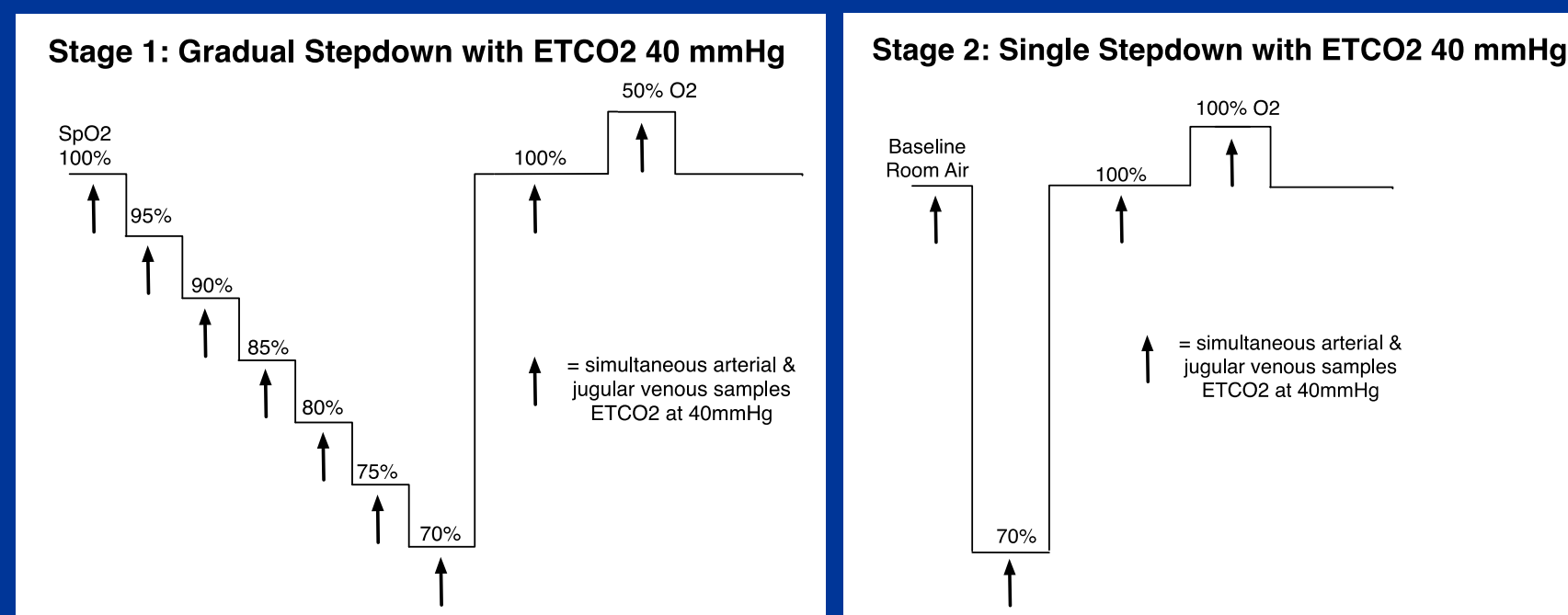

Each step was maintained for 6 minutes to establish steady-state end-tidal CO<sub>2</sub> and O<sub>2</sub> tensions during which a single jugular bulb blood sample and two arterial blood samples were drawn simultaneously.

Co-oximetry was used to determine the jugular venous saturation (SjvO<sub>2</sub>) and arterial saturation (SaO<sub>2</sub>). The two arterial samples per plateau were averaged to provide one SaO<sub>2</sub> value per plateau.


The arteriovenous (SavO<sub>2</sub>) saturation was then calculated as a 70:30 ratio of SjvO<sub>2</sub> to SaO<sub>2</sub>.

$$SavO_2 = [0.70 \times SjvO_2] + [0.30 \times SaO_2]$$

A predetermined ARMS value of 5% was set as the criterion for acceptable accuracy.

## Results

24 subjects completed the study (calibration=13; validation=11). Readings were obtained in all subjects and no subjects were excluded from the analysis. Absolute accuracy of rSO<sub>2</sub> compared to calculated SavO<sub>2</sub> as measured by ARMS was 4.1%.

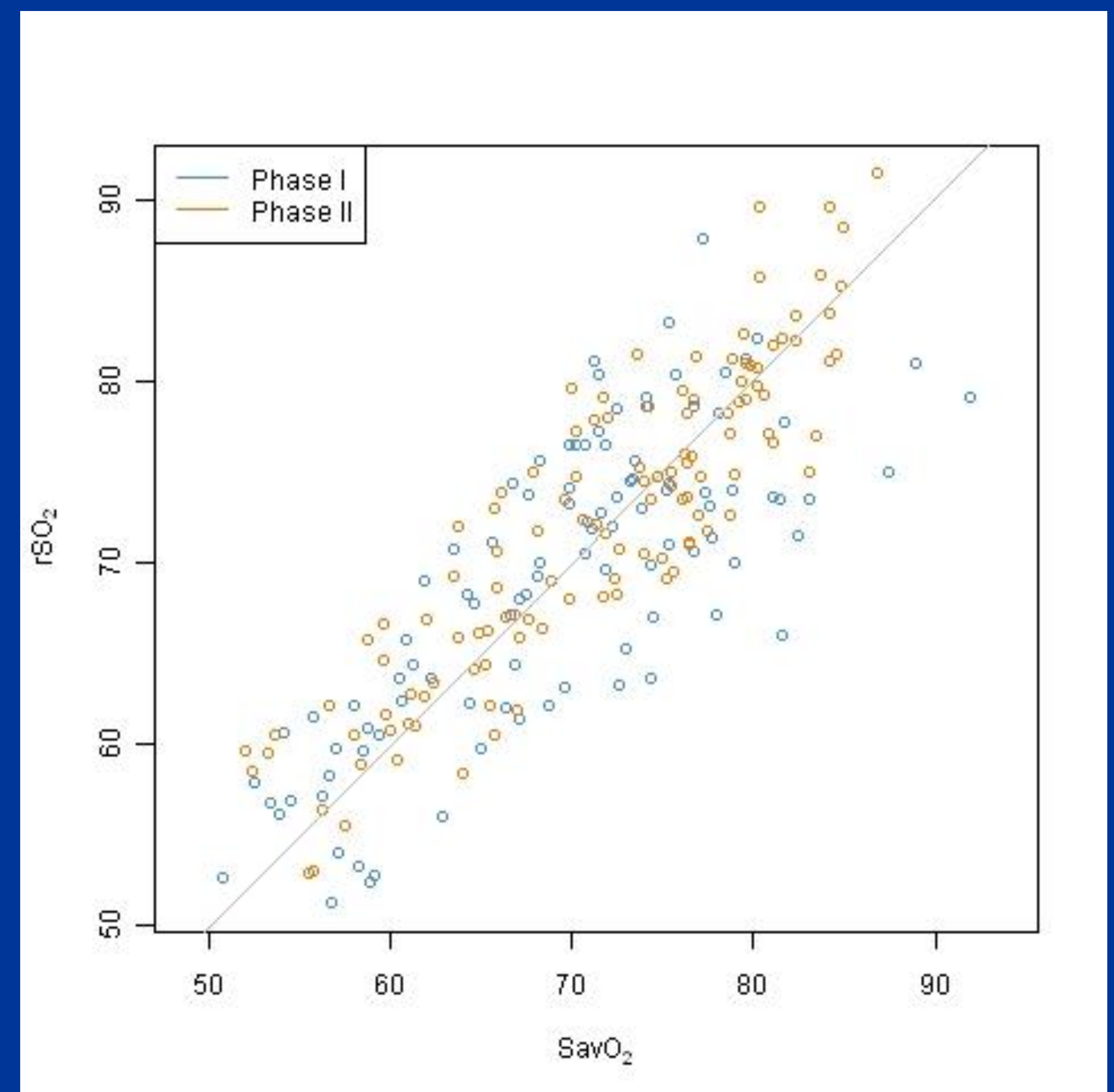

## Conclusion

This study confirms that Nonin's 4-wavelength cerebral oximeter provides an accurate measure of the calculated cerebral tissue oxygen saturation during deliberate oxygen desaturation in healthy volunteers.

DUKE

Extraordinary Care – Through a Culture of Innovation